
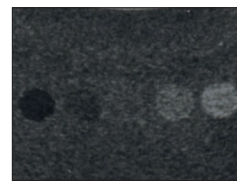



Ultrasound QA Phantom (Multi Purpose Phantom N-365)

Durable and stable

Useful both for daily assessment and further research. Gray scale for evaluation, cyst targets with non-resonance cylinders, line targets for geometrical evaluation including close range (dead zone) resolution, axial and lateral resolution are prepared for scanning. The phantom is designed to allow scanning from all four side walls.

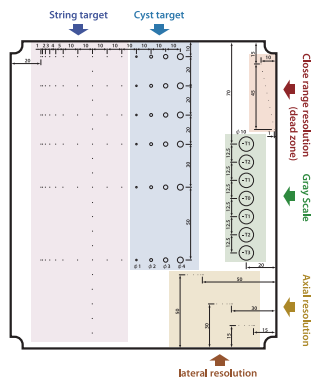

Gray Scale

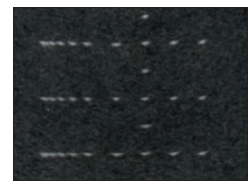

String target

SPECIFICATIONS

Size: 19 x 22 x 7 cm
7.6 x 8.8 x 2.8 in.

Sonic velocity: 1432m/ sec at 25 degrees C
Attenuation coefficient: 0.57 dB/ cm MHz at 25 degree C
Acoustic impedance: 1.38 rayl at 25 degree C

SET INCLUDES

- 1 phantom
- 1 carrying case

DESCRIPTIONS

MATERIALS

Urethane elastomer, acryl, nylon
Latex free

REPLACEMENT PARTS

RECOMMENDED DEVICES

-