


# Pro-MAM Geometry

03-107

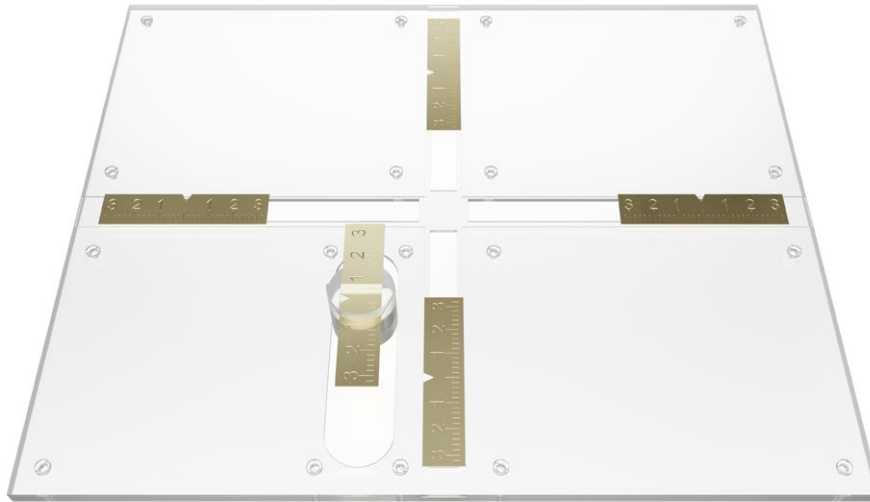

This test device was developed to provide the user with a convenient, quick and accurate way of testing collimation of mammography units. It is very simple to use. Instead of using all those hard-to-find coins, you only use our test tool.

## Technical data (can be modified to customer specifications):

- dimensions: 300 x 240 mm
- measurement can be quickly and easily repeated
- compression paddle rests on peg exactly 4.2 cm above the Bucky
  - no measurement of compression paddle height needed
- adaptable for 18 x 24 cm, 24 x 30 cm and magnification stand testing

## Product features:

- Complies with:
  - IEC 61223-3-2
  - MQSA testing requirements as contained in the ACR Mammography QC Manual
- the Manual provides detailed guidelines for carrying out each test, results assessment and registration